




Grand Master BTE 13 Information Sheet

	BTE 13	BTE 13
	Power	Power (Dir)
		
BTE Product Options		
Battery Size	13	13
Directional	✗	✓
Multi Program	✓	✓
Multi Function Digital VC	✓	✓
Telecoil	✓	✓

Fitting Range

Grand Master BTE 13 - Power

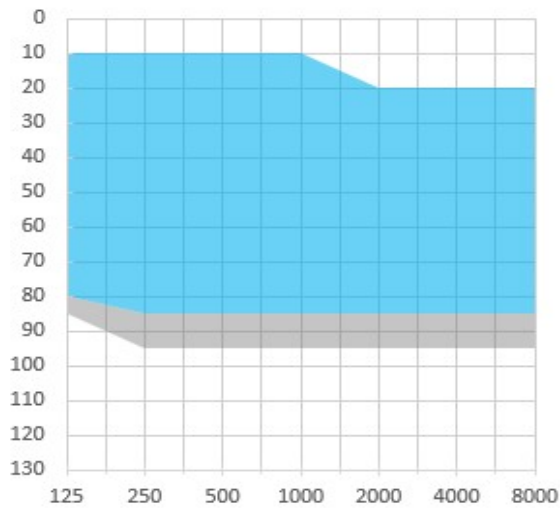


Technical Description


1. Ear Hook
2. Multi-Function Digital VC & Program Change
3. Battery Compartment
4. Mic(s)
5. Option door



Programming via cable 4 pin

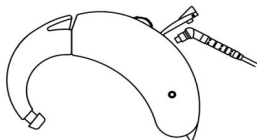

1. Open the option door for access to the programming socket.

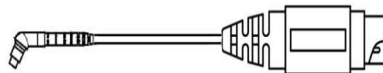

2. Connect the programming plug into the socket.
Use the programming cable 4 Pin only.

3. Please ensure that the hearing instrument is on during programming.