

Grand Master

Product Information

version 1.0





Grand Master Features Overview

	BTE 675	BTE 675	BTE 13	BTE 13
	SP (Dir)	SP	Power (Dir)	Power
12 Gain Bands	✓	✓	✓	✓
6 WDRC Channels	✓	✓	✓	✓
4 Listening Programs	✓	✓	✓	✓
Ability to disable additional program(s)	✓	✓	✓	✓
Accurate Fitting with Multilingual ALPSiFit	✓	✓	✓	✓
Adaptive Feedback Cancellation	✓	✓	✓	✓
Directionality	Adaptive Directionality	Omnidirectional	Adaptive Directionality	Omnidirectional
Adaptive Noise Reduction	✓	✓	✓	✓
Auto save last used program & Volume Control position	✓	✓	✓	✓
Clarifier-Multichannel	✓	✓	✓	✓
High Fidelity Sound	✓	✓	✓	✓
Hearing Aid Ear Indicator (Blue for Left & Red for Right)	✓	✓	✓	✓
Low Battery Alert (Customizable Tone Frequency & Loudness)	✓	✓	✓	✓
Low Battery Consumption	✓	✓	✓	✓
Multi Function Digital Volume Control	✓	✓	✓	✓
Program Change Alert (Customizable Tone Frequency & Loudness)	✓	✓	✓	✓
Program switch & Volume Control deactivation	✓	✓	✓	✓
Special Setting for Covid-19 Masks	✓	✓	✓	✓
Startup Settings - Customizable	✓	✓	✓	✓
Telecoil facility (Telecoil Gain – Programmable)	✓	✓	✓	✓
VC Range Customizable	✓	✓	✓	✓
Water repellent coating Available	✓	✓	✓	✓
Wireless Connectivity Bluetooth & FM (optional)	✓	✓	✓	✓




Ashirwad Complex, D-1 Green Park, New Delhi - 110016 (India) Email: sales@alps.in Web: www.alps.in

Grand Master BTE 13 Key Features

	BTE 13	BTE 13
	Power	Power (Dir)
		
BTE Product Options		
Battery Size	13	13
Directional	✗	✓
Multi Program	✓	✓
Multi Function Digital VC	✓	✓
Telecoil	✓	✓

Fitting Range

Grand Master BTE 13 - Power

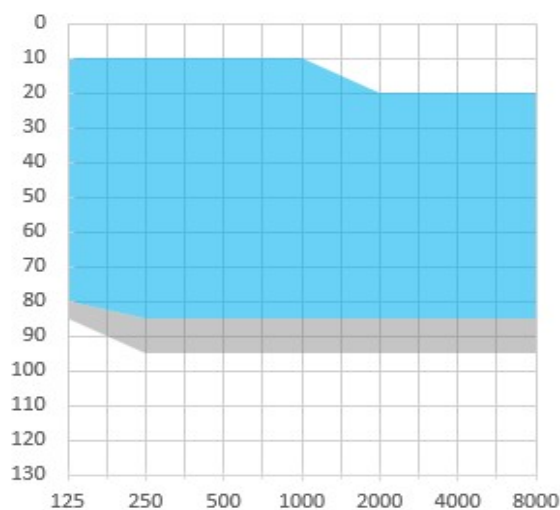


Technical Description


1. Ear Hook
2. Multi-Function Digital VC & Program Change
3. Battery Compartment
4. Mic(s)
5. Option door



Programming via cable 4 pin

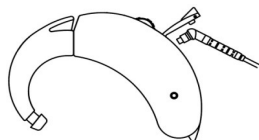

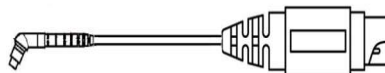

1. Open the option door for access to the programming socket.



2. Connect the programming plug into the socket.
Use the programming cable 4 Pin only.


3. Please ensure that the hearing instrument is on during programming.

Grand Master BTE 675 Key Features

	BTE 675	BTE 675
	SP	SP (Dir)
		
BTE Product Options		
Battery Size	675	675
Directional	✗	✓
Multi Program	✓	✓
Multi Function Digital VC	✓	✓
Telecoil	✓	✓


Fitting Range

Grand Master BTE 675 (SP)

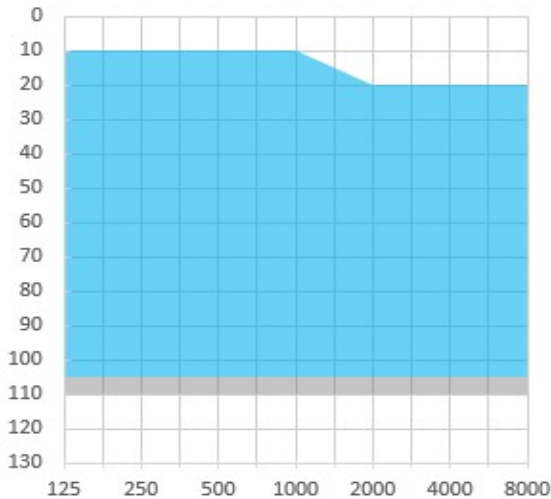


Technical Description


1. Ear Hook
2. Multi-Function Digital VC & Program Change
3. Battery Compartment
4. Mic(s)
5. Option door

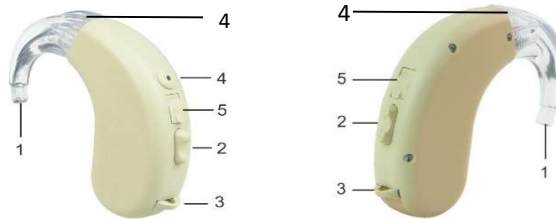

Programming via cable 4 pin

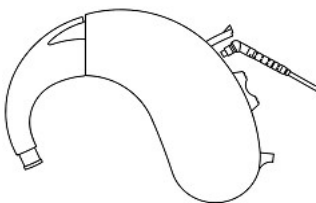

1. Open the option door for access to the programming socket.

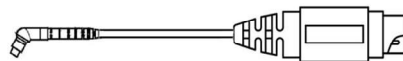

2. Connect the programming plug into the socket.
Use the programming cable 4 Pin only.
3. Please ensure that the hearing instrument is on during programming.